




## INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification n 5 : <b>C12Q 1/04, 1/34</b>		A1	(11) International Publication Number: <b>WO 90/12888</b>
			(43) International Publication Date: <b>1 November 1990 (01.11.90)</b>
(21) International Application Number: <b>PCT/US90/02346</b>		(74) Agents: MAKI, David, J. et al.; Seed and Berry, 6300 Columbia Center, Seattle, WA 98104-7092 (US).	
(22) International Filing Date: <b>27 April 1990 (27.04.90)</b>			
(30) Priority data: <b>345,033 27 April 1989 (27.04.89) US</b>		(81) Designated States: AT (European patent), AU, BE (European patent), CA, CH (European patent), DE (European patent), DK (European patent), ES (European patent), FR (European patent), GB (European patent), IT (European patent), JP, LU (European patent), NL (European patent), SE (European patent).	
(71) Applicant: BIOCONTROL SYSTEMS, INCORPORATED [US/US]; 19805 North Creek Parkway, Bothell, WA 98011 (US).			
(72) Inventors: WARD, Robert, N., Jr. ; 4237 Second Avenue N.E., Seattle, WA 98105 (US). DESROSIER, John, P. ; 3520 Burke Avenue N., Seattle, WA 98103 (US). MARSHALL, Elliott, D., III ; 1411 Grant Avenue S., Apt. #M-305, Renton, WA 98055 (US). FORD, Judith ; 17939 95th Place N.E., Apt. #6, Bothell, WA 98011 (US). MALLINAK, Nancy, J., S. ; 1900 N.W. 89th, Apt. #2, Seattle, WA 98117 (US).		Published <i>With international search report. Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i>	

## (54) Title: PRECIPITATE TEST FOR MICROORGANISMS

## (57) Abstract

An improved method for testing for the presence of a particular microorganism and/or group of microorganisms characterized by a particular enzyme is provided. The method uses a dye-forming and/or fluorogenic-generating substrate throughout a polymer matrix or on the surface of a solid support member that forms a colored precipitate or fluorescent by-product when cleaved by the enzyme. The precipitate is concentrated on the polymer matrix and/or solid support member to create a visible reaction product, whereas the fluorescent by-product will normally be dispersed into the surrounding medium to some extent. The amount of dye-forming substrate needed is independent of the sample size. The present invention further comprises an enzyme indicator device comprising a dye-forming and/or fluorogenic-generating substrate throughout a polymer system or on the surface of a solid support member.

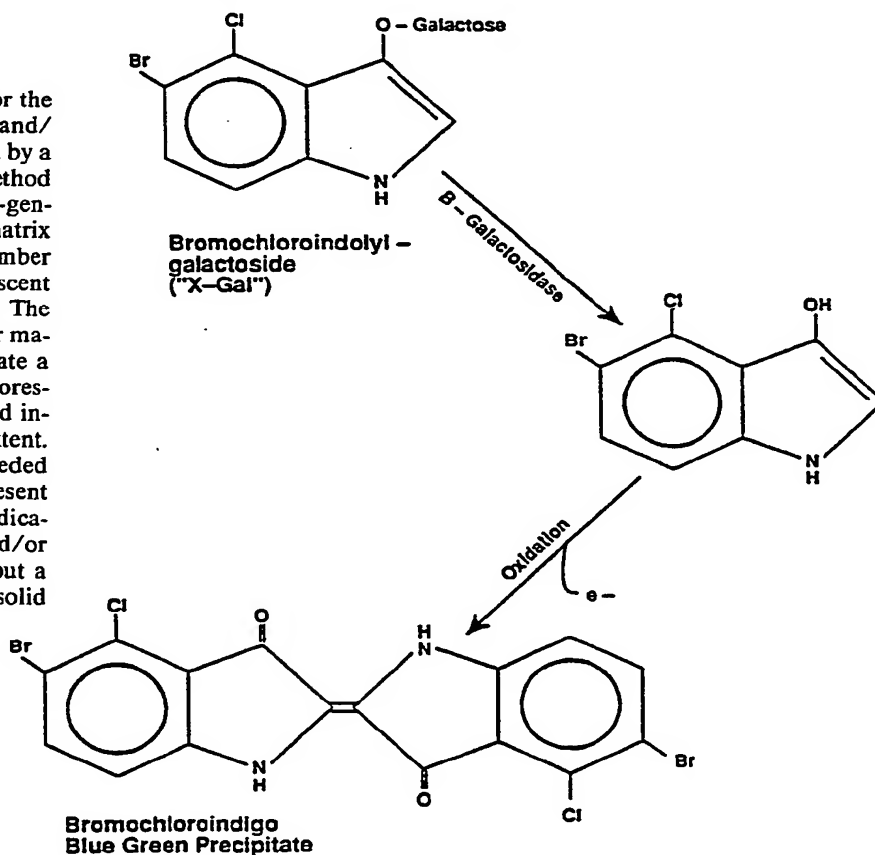

### DESIGNATIONS OF "DE"

Until further notice, any designation of "DE" in any international application whose international filing date is prior to October 3, 1990, shall have effect in the territory of the Federal Republic of Germany with the exception of the territory of the former German Democratic Republic.

#### FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	ES	Spain	MC	Monaco
AU	Australia	FI	Finland	MG	Madagascar
BB	Barbados	FR	France	ML	Mali
BE	Belgium	GA	Gabon	MR	Mauritania
BF	Burkina Faso	GB	United Kingdom	MW	Malawi
BG	Bulgaria	GR	Greece	NL	Netherlands
BJ	Benin	HU	Hungary	NO	Norway
BR	Brazil	IT	Italy	RO	Romania
CA	Canada	JP	Japan	SD	Sudan
CF	Central African Republic	KP	Democratic People's Republic of Korea	SE	Sweden
CG	Congo	KR	Republic of Korea	SN	Senegal
CH	Switzerland	LI	Liechtenstein	SU	Soviet Union
CM	Cameroon	LK	Sri Lanka	TD	Chad
DE	Germany, Federal Republic of	LU	Luxembourg	TG	Togo
DK	Denmark			US	United States of America

Description

## PRECIPITATE TEST FOR MICROORGANISMS

5

Technical Field of the Invention

The present invention relates to an improved method for assaying for a particular microorganism or group of microorganisms wherein the particular microorganism or group of microorganisms are characterized by the presence of a particular enzyme. The assay uses a dye-forming substrate within a polymer matrix or on the surface of a solid support. The dye precipitates when cleaved by the enzyme. The present invention further comprises an enzyme indicator device comprising a dye-forming substrate, polymer system, and a solid support.

15 Background of the Invention

Human gastrointestinal disease can be caused by a variety of microorganisms. Common vehicles for infection are contaminated foods and water. To reduce the incidence of such disease, foods and water destined for human consumption are routinely tested for their sanitary quality. Instead of testing for a multitude of different enteric pathogens, laboratories test for the presence of indicator organisms. Coliform bacteria are typically used as the primary indicators of sanitary quality because they are commonly associated with the gastrointestinal tracts of warm-blooded animals. The presence of coliform bacteria, especially Escherichia coli, in high numbers in foods or water suggests that there may be fecal contamination, and contraindicates human consumption.

Coliform bacteria are distinguished from other organisms and from their close relatives in the family Enterobacteriaceae by their ability to ferment lactose to acidic and gaseous (CO<sub>2</sub> and H<sub>2</sub>) end products. Certain non-coliform bacteria may ferment lactose to the same end products, but the growth of these organisms in the coliform assay is usually minimized or prevented by the selective properties of the bacteriological media used for testing.

Typically, the presence of coliform bacteria in a sample is determined by adding the material to a liquid bacteriological growth medium and incubating the mixture at a temperature which is conducive to bacterial growth. Incubation of the mixture of sample and growth medium is important because the assay must detect low levels of coliform bacterial contamination in the sample. Incubation of the mixture results in the multiplication of coliforms to a level of

approximately  $10^6$  to  $10^8$  cells per milliliter of culture, wherein their presence may be detected by any of a number of techniques. The variety of different liquid bacteriological growth media which have been used for coliform detection share two common properties: ~~they contain the disaccharide lactose and they also~~

- 5 contain chemical agents which selectively inhibit the growth of non-coliform microorganisms. Selection is important because the sample invariably contains a variety of microorganisms, and the success of the assay depends on the coliforms not being overgrown by the non-coliforms.

Coliform assays can be performed using either solid or liquid media.

- 10 Assays which employ a solid medium allow a viable cell count. Sample is added to the solid medium and discrete colonies of coliforms are enumerated. Alternatively, samples may be added to liquid media. Coliforms are detected through the formation of characteristic metabolic end products. The liquid medium format may be qualitative or quantitative. The liquid medium format is  
15 preferred for samples containing fewer coliforms (e.g., less than 10 organisms per milliliter), or samples containing particulate material (e.g., food or dairy samples) which obscures colony visualization.

- Most coliform assays in a broth or agar take place in two discrete stages: presumptive and confirmed. First, a presumptive assay provides an  
20 indication of possible coliform presence. In the confirmed stage, presumptively positive cultures or typical colonies are subcultured into a second, more selective medium. In principle, the confirmed medium eliminates false positive results. Together, the two stages of the assay require 48 to 96 hours for completion. Therefore, there is a need in the art to provide accurate results in a more timely  
25 manner.

- Generally, detection is based upon the end products of metabolic pathways of coliform bacteria. For example, acid production by coliform bacteria is generally detected through the incorporation of pH indicators in the medium. Acid formation can be detected by a change in color from purple to yellow of the  
30 indicator bromocresol purple in a liquid medium such as Clarke's medium. Acid production is detected by the formation of colonies which are dark-centered due to the reaction between the acid and the indicator neutral red in a solid medium such as violet red bile agar (VRBA) and MacConkey agar.

- Further, gaseous end products are usually detected by the presence  
35 of gas bubbles in a liquid. Gas bubbles may be entrapped in a smaller, inverted test tube or an inverted vial within the culture tube or in a special portion of the culture device. For example, the BioControl ColiTrak™ product entraps gas

bubbles in a dome associated with the top portion of the device. Petrifilm™ (3M), with a solid growth medium, entraps gas bubbles in close proximity to a bacterial colony. Alternatively, gaseous end products may be detected by nonvisual means such as by electrochemical detection of hydrogen, by radiometric detection of  $^{14}\text{CO}_2$  released from the fermentation of radiolabelled lactose, or by impedimetric or gas chromatographic detection of organic compounds produced during fermentation.

An alternative approach to coliform detection is based upon the detection of coliform-associated enzyme activity rather than metabolic end products. Generally, an enzyme assay approach can yield quantitative estimates for coliform bacteria comparable to confirmatory results using lengthy culture testing procedures. The enzyme beta-galactosidase is a bacterial enzyme used for the fermentation of lactose. Beta-galactosidase hydrolyzes lactose to its component sugars glucose and galactose. Coliform bacteria typically express this enzyme. For example, Warren et al., Appl. Environ. Microbiol. 35:136-41, 1978, refers to a coliform testing method which uses the chromogenic beta-galactosidase substrate, o-nitrophenyl-beta-D-galactoside (ONPG), for quantitating fecal coliforms in water. In Warren et al., the time of appearance of the yellow reaction product o-nitrophenol (ONP) is inversely related to the initial concentration of coliforms in the test sample.

The product Colilert™ (Access Medical) uses the ONPG enzyme indicator for detecting coliforms in drinking water. The water sample is used to solubilize a basal growth medium containing ONPG. Beta-galactosidase from the bacteria growing in the culture hydrolyzes the ONPG to produce galactose, which serves as the sole carbon and energy source for growth. ONP formation indicates enzyme activity. After 24 hours of incubation, the presence of a yellow color throughout the liquid culture of 10 milliliters signifies that coliform bacteria were present in the water sample.

The fluorogenic beta-galactosidase substrate fluorescein-beta-D-galactoside is another indicator of activity. Cundell et al., Proc. Water Reuse Symp. 3:1895-99, 1979, refers to an assay incubating coliform bacteria in a medium containing the fluorescein-beta-D-galactoside. Coliform bacteria are determined quantitatively by flow cytometry. The fluorescein moiety of the substrate becomes concentrated within the cell after cleavage by beta-galactosidase and imparts fluorescence to the cells under ultraviolet illumination. U.S. Patent No. 4,242,447 refers to a fluorescent assay for enumerating coliforms using a fluorogenic substrate. Specifically, U.S. Patent No. 4,242,447 refers to a process emulsifying a

water sample and a fluorogenic substrate with an oil to form oil droplets containing coliform bacteria. Beta-galactosidase activity cleaves the fluorogenic substrate and forms a fluorescent oil droplet which is counted in a fluorescence detector. Cundell et al., Dev. Ind. Microbiol. 20:571-77, 1979, refers to a similar technique that yields oil-encapsulated cells on a microscope slide. Fluorescent droplets are counted under a fluorescence microscope. The Cundell method has been used with sewage samples. The sensitivity of the Cundell method is  $10^5$  cells per milliliter.

The bacterium, E. coli, is a coliform which may be assayed separate and apart from the coliform bacteria. E. coli presence is considered to be a more reliable indicator of fecal contamination. Additionally, certain strains of E. coli are pathogenic for humans and animals. The standard E. coli detection method subcultures positive presumptive cultures into EC broth and incubates at 45.5°C. Gas-positive EC cultures are then streaked onto a differential agar medium. The isolates are identified by biochemical characterization.

E. coli is a relatively unique organism by possessing the enzyme beta-glucuronidase. Consequently, the fluorogenic substrate 4-methylumbelliferyl- $\beta$ -D-glucuronide (MUG) is used to detect E. coli in food and water samples. The E. coli detection procedure consists of inoculating cultures and incubating by standard methods in the presence of MUG. Fluorescence which develops during 24 hours incubation indicates the presence of E. coli because beta-glucuronidase from E. coli cleaves MUG to a fluorescent product. Thus, the use of MUG-containing media can shorten the detection time for E. coli. Unfortunately, MUG is an expensive reagent and large sample sizes, such as water samples, require large quantities of MUG and make the test prohibitively expensive. Currently, the recommended concentration of MUG is 50-100 micrograms per milliliter of final culture. Typical volumes for food and water tests are 90 and 100 milliliters respectively. Thus, the cost of the MUG reagent makes it unattractive for food testing and for use in the PA type of water testing.

One can assay for beta-galactosidase and beta-glucuronidase activities as a means for determining the presence of coliforms and E. coli, respectively. Enzyme detection assays offer advantages over the detection of fermentation end products. First, the enzyme assay approach is more sensitive. One enzyme can cleave many substrate molecules. Each molecule of substrate cleaved by an enzyme yields a fluorescent or colored reporter product. Therefore,

the signal is amplified by an enzyme. With certain substrates, as few as  $10^4$  metabolically active cells can yield a signal within 24 hours.

By contrast, it takes approximately  $10^7$  to  $10^8$  metabolically active cells to produce a visible gas bubble within 24 hours, because beta-galactosidase cleavage of lactose does not result in stoichiometric amounts of acid and gas end products due to physiological regulation of metabolism. The majority of the carbon in the lactose molecule is incorporated into cellular biomass and is not used to regenerate oxidized NAD during fermentation. The bacterial fermentation enzymes responsible for acid and gas production are at the end of the fermentation pathway and are tightly regulated. For instance, pyruvate-formate lyase is inactive under aerobic conditions and is activated only when the culture becomes anaerobic. Formate dehydrogenase, a key enzyme in gas production, is not synthesized in the presence of oxygen or when the culture pH is 6 or greater. Hydrogenase activity and synthesis are negatively regulated by alternative oxidants such as oxygen, nitrate and nitrite. Accordingly, it is possible that growth conditions of anaerobiosis, acidic pH and the absence of external oxidants are either not achieved or are achieved only slowly. Moreover, much of the lactose in the culture will have been cleaved by the time the acid- and gas-producing enzymes become active. Therefore, enzyme assays for beta-galactosidase circumvent the fermentation pathway regulation problems and yield a reporter group with every substrate cleavage event.

The presence of anaerogenic strains of genera which belong to the coliform group further complicates gas detection assay procedures. These strains possess beta-galactosidase activity, but do not produce gas from lactose. Thus, only an assay for beta-galactosidase activity would detect the anaerogenic strains of the coliform group.

It is possible to detect coliform-associated enzymes within 24 hours of incubation. Beta-galactosidase and beta-glucuronidase assay results correlate with the confirmed presence of coliforms and E. coli. When aliquots from 24 hour enzyme assay-positive coliform cultures are streaked onto a differential medium, such as eosin methylene blue agar, and the agar incubated, typical coliform and E. coli colonies are almost always recovered. Members of the coliform genera are the taxonomically identified colonies from beta-galactosidase-positive cultures. E. coli is the taxanomically identified organism from beta-glucuronidase-positive cultures. These observations indicate that the first incubation (presumptive) stage of the standard coliform assay can be shortened from 48 to 24 hours, and that the second (confirmed) incubation stage is unnecessary. Consequently, enzyme assays

for coliforms and E. coli can reduce the time required to produce a confirmatory result by as much as 72 hours, and simultaneously, provide significantly better sensitivity. However, the cost is more expensive due to the reagents, especially for large sample sizes.

5 Food and dairy samples present problems detecting coliforms. There is a large diversity of foods and the physical/chemical nature of some makes the detection of either acidic or gaseous end products or ~~certain types of enzyme assays difficult.~~ For example, gas bubbles may be difficult to detect in a  
10 broth due to the turbidity from a food sample such as nonfat dried milk. The color associated with a food product such as dried cheese, chili powder or green or yellow vegetables can obscure the color of a pH indicator or a colored reporter, such as ONP, in a beta-galactosidase assay using ONPG as a substrate. Further, a food may be inherently acidic, or it may internally trap gas bubbles during sample preparation. Finally, food chemistry may interfere with an assay. Foods add  
15 additional nutrients to a culture broth and these nutrients may be metabolized by non-coliform bacteria to produce gaseous or acidic end products leading to a false-positive result. For example, cake mixes or other sweetened products contain significant amounts of disaccharides and monosaccharides, which can be converted to acid and/or gas by non-coliforms.

20 It is important to choose an appropriate substrate for beta-galactosidase and beta-glucuronidase assays for coliform bacterial and E. coli detection. For example, beta-galactosidase cleaves the substrate ONPG to a yellow product, ONP. However, certain colors inherent in or added to food products may obscure the yellow color of the ONP reporter product. Food  
25 substances which can obscure the yellow color include hemoglobin in red meats, chlorophyll in leafy greens and vegetables, beta-carotene and other yellow and orange pigments in fruits and vegetables, natural pigments in spices, and artificial and natural colorants. Turbid water and water with a high humic content can also obscure weakly yellow positive ONP reactions. Moreover, standard  
30 bacteriological growth media ingredients, such as protein or yeast hydrolysates, impart a yellow or golden color to the media and can obscure the presence of ONP. For example, the standard method coliform medium, lauryl sulfate tryptose (LST) broth, has a golden hue. The result of performing an ONPG-type enzyme assay in LST broth and other supplemented media will lead to a significant  
35 number of false-negative results.

An additional problem associated with the use of ONPG is the hazardous nature of the ONP cleavage product. ONP is a volatile compound and



a known respiratory eye and skin irritant. The presence of ONP makes activities associated with handling of cultures and washing of glassware hazardous to lab workers. There are also toxic waste disposal problems when using ONP.

Similarly, beta-glucuronidase assays use MUG, which generates a  
5 fluorescent signal. Pigments associated with foods have been observed to quench the fluorescent signal. For example, the pigment in chili powder obscures the fluorescence reaction of the reporter product. A further problem with MUG is that the pH range for maximal fluorescence is outside of the pH range produced during lactose fermentation. MUG is maximally fluorescent above pH 10. Yet,  
10 acids produced from lactose fermentation typically lower the coliform culture medium to pH 5 to 6, obscuring fluorescence reactions. Finally, as noted above, beta-glucuronidase assays which employ chromogenic substrates may be subject to interference from some types of samples. Therefore, current enzyme assays for beta-glucuronidase have pH and interference problems from the food samples,  
15 and are inappropriate for use in certain types of foods.

Fluorogenic substrates, such as MUG have been used in agar to detect E. coli in food samples. The fluorescent cleavage product, 4-methylumbelliferone, is soluble and will diffuse from E. coli colonies in agar (a solid medium) to add to the difficulties associated with the determination of which  
20 colonies contain the enzyme beta-glucuronidase and which do not.

An additional complication is that the final volume of culture broth for different types of samples varies from 1 ml to 1 liter. For each test, the molarity of the chromogenic or fluorogenic substrate within the culture medium would need to be constant. This would require consumption of large amounts of  
25 expensive test reagents. It would be advantageous to a food processor if all coliform testing could be done in a single assay format which features easily read results across all test samples, regardless of size, and which economizes on test reagents.

A typical food processing company or water treatment facility tests  
30 incoming and outgoing foods, process water, waste water or environmental samples. Therefore, an ONPG- or MUG-based assay may only be appropriate for certain types of samples. Thus, there is a need in the art for assays which are functional for a variety of samples, such as foods, turbid and/or colored water and environmental samples and which eliminate the difficulties associated with  
35 correlating different test methodologies.

The indigogenic substrate 5-bromo-4-chloro-3-indolyl-beta-D-galactoside (X-Gal) has been used in solid media for the detection of E. coli

colonies which express beta-galactosidase activity. Beta-galactoside cleaves X-Gal to galactose and the indolyl derivative. The indolyl derivative dimerizes to form a substituted indigo, which has an intense blue color.

- Frampton et al., J. Food Prot. 51:402-04, 1988, refers to a peptone-tergitol agar supplemented with 5-bromo-4-chloro-3-indolyl- $\beta$ -D-glucuronide (X-Gluc) to differentiate E. coli from other bacterial colonies in artificially inoculated raw minced chicken. Ley et al., Ann. Meet. Am. Soc. Microbiol., Abstract Q35:288, 1988, refers to a method to enumerate E. coli in water and sewage using a membrane filtration technique with a glycerol medium and 3-indolyl-beta-D-glucuronide substrate. E. coli colonies appeared blue in the semisolid medium. Watkins et al., Appl. Environ. Microbiol. 54:1874-75, 1988, refers to a pour plate method using X-Gluc for detecting E. coli in shellfish and waste water.

A product Petrifilm™ (3M) uses a solid medium and X-Gluc for detecting E. coli. E. coli colonies appear blue.

- Because of the high cost of the indogenic compounds X-Gal and X-Gluc, for the detection of beta-galactosidase and beta-glucuronidase activity have not found widespread use or have been used with relatively small culture volumes. Therefore, there is a need in the art for a beta-galactosidase test system that can economically use small quantities of X-Gal and X-Gluc for food and water testing, because the intense blue color will solve many of the problems associated with the yellow color of ONP and other colored reporters.

- In summary, there is a need in the art, therefore, for a coliform detection assay which offers the improved sensitivity and speed associated with enzyme assays, the capacity to detect anaerogenic coliforms, successful applications to food, water and environmental samples, safety of handling and disposal, and economy in the consumption of expensive test reagents.


### Summary of the Invention

- The present invention comprises a method for assaying an enzyme associated with a particular microorganism or a group of microorganisms and an enzyme indicator device that is useful for assaying for a particular microorganism and/or group of microorganisms. In one aspect, the method comprises adding sample and an enzyme indicator device to growth medium, incubating at a permissive temperature and reading the results on the enzyme indicator device. Preferably, the growth medium is a liquid. Another method comprises adding a sample to a growth medium; incubating the sample and growth medium at a permissive temperature for the growth of microorganisms until there is evidence

of microorganism growth within the growth medium; adding an enzyme indicator device to the growth medium, wherein the enzyme indicator device comprises a dye-forming and/or fluorogenic-generating substrate in a polymer system on a solid support, or a dye-forming or fluorogenic-generating substrate dried on a solid support, and wherein the dye-forming substrate is cleaved by the enzyme to form a colored precipitate or fluorescent reaction; incubating the sample, growth medium, and enzyme indicator device at a permissive temperature and in an oxidizing environment from about one hour to about forty-eight hours; and visually detecting the presence of the particular microorganism or group of microorganisms by a color change or fluorescent reaction on the enzyme indicator device. Preferably, an oxidizing agent is added to the growth medium along with the enzyme indicator device. The enzyme indicator device preferably comprises a dye-forming and/or fluorogenic-generating substrate in a polymer system. Alternatively, the enzyme indicator device comprises a dye-forming and/or fluorogenic-generating substrate adhering to the surface of a solid support. More particularly, the enzyme indicator device further comprises a solid support member, wherein the dye-forming and/or fluorogenic-generating substrate in the polymer system is layered onto a surface of the solid support system. It is preferred that the solid support floats on the surface of a liquid so as to have an oxidizing environment from air without the need for additional oxidizing agents.

Within another aspect of the present invention, a method for assaying for an enzyme associated with a particular microorganism or a group of microorganisms is provided, comprising adding a sample to a growth medium; adding a dye-forming substrate in a polymer system, wherein the dye-forming substrate forms a precipitate in association with the polymer, when cleaved by the enzyme; contacting the dye-forming substrate in the polymer with the sample and the growth medium; incubating the sample, growth medium and dye-forming substrate in polymer at a permissive temperature for the growth of particular microorganism or group of microorganisms; and detecting the presence of the particular microorganism or group of microorganisms by a colored precipitate in association with the polymer.

The dye-forming substrate forms an insoluble precipitate upon enzymatic cleavage. Preferably, the dye-forming substrate is a substituted indigo derivative of the formula:

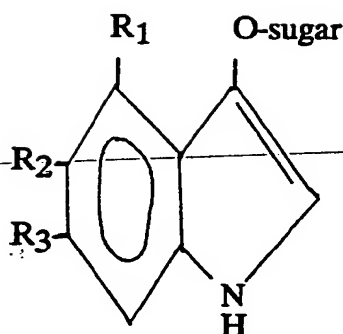

10

wherein each R substituent is a halo group or hydrogen. The halo group is a chloro, bromo, or iodo group. When R<sub>1</sub>, R<sub>2</sub> and R<sub>3</sub> are H, the precipitate dye formed is indigo.

15

It is preferred that the dye-forming substrate is insoluble in aqueous systems so as to prevent the dye-forming substrate from leaching out of the polymer or from coming off the surface of the solid support when exposed to the growth medium. In this way, it is possible to conserve costs by reducing the amount of dye-forming substrate needed for each test. Thus, the amount of dye-forming substrate used is independent of the test size or the sample size.

20

Moreover, the analysis of particular microorganisms or group of microorganisms by enzyme assays provides a more rapid and more sensitive assay system than conventional growth and pure culture detection methods.

The advantages of the inventive method and the use of the inventive enzyme indicator device allow the dye-forming substrate to be deployed on a defined reaction surface rather than dispersed throughout a solution. The reaction surface is the surface of a solid support member or a polymer matrix. This results in an assay format which is independent of culture volume, thus leading to an economy of reagent consumption. Second, the insoluble indigo dye-forming substrates from the indigogenic substrates range in color from blue to purple, colors which are not masked or quenched by pigments or turbidity associated with foods or certain types of water. Thus, the indigogenic substrates are amenable to detection of coliform bacteria in both food and water samples. Third, the indigogenic substrates offer the speed and sensitivity associated with enzyme assays. Fourth, the assay for beta-galactosidase eliminates the problems associated with enzyme-regulation in the lower portion of the lactose fermentation pathway. Fifth, an enzyme assay ensures that signal is being generated in the first step of the pathway every time the enzyme performs a hydrolysis of the dye-

35

forming substrate, rather than at the end of the pathway where only a portion of the cleaved galactoside ends up.

Within another aspect of the present invention, the method for assaying for an enzyme associated with a particular microorganism and a group of microorganisms comprises (a) adding a sample to a liquid growth medium; (b) adding a dry-forming or fluorogenic-generating substrate on a defined reaction surface wherein the substrate forms a precipitate or a fluorescent by-product when cleaved by the enzyme; (c) incubating the sample, growth medium and substrate at a permissive temperature for the growth of the particular microorganism and group of microorganisms; and (d) detecting the presence of the particular microorganism and group of microorganisms by a colored precipitate or fluorescence on the defined reaction surface or in the surrounding medium. Within a preferred embodiment, the fluorogenic-generating substrate is 4-methylumbelliferyl- $\beta$ -D-glucuronide.

#### Brief Description of the Drawing

The Figure shows the reaction pathway for bromochloroindolylgalactoside (X-Gal) hydrolysis to form bromochloroindigo. Bromochloroindolylgalactoside is first cleaved by beta-galactosidase from the particular microorganism or group of microorganisms to remove the galactose group and form bromochloroindoxyl, which then dimerizes by an oxidation reaction to the highly colored blue-green precipitate, bromochloroindigo. Oxidation can be achieved either by an oxidizing agent in solution or by exposure of bromochloroindoxyl to air.


#### Detailed Description of the Invention

The present invention is a method for assaying an enzyme associated with a particular microorganism and/or a group of microorganisms. Preferably, the particular microorganism is E. coli and the group of microorganisms is the coliform bacteria. The enzyme beta-galactosidase is characteristic of the coliform bacteria and beta-glucuronidase activity is characteristic of E. coli. The dye-forming substrates preferably are mixed in a polymer system and form an insoluble precipitate upon enzymatic cleavage. Alternatively, the dye-forming substrates may adhere to a solid support member which also is the reaction surface for forming the insoluble precipitate. Preferably, the dye-forming substrate is insoluble in an aqueous liquid, such as a growth medium so as to avoid dispersion of the dye-forming substrate throughout

the growth medium and instead, have the dye-forming substrate concentrated within the polymer system or on the surface of the solid support member.

The reaction of the dye-forming substrate to form a highly colored precipitate requires the presence of the specific enzyme to cleave off the sugar group and an oxidizing environment to form the highly colored dimer precipitate from the cleaved dye-forming substrate. The oxidizing environment can be created by having the dye-forming substrate in the polymer exposed to air or by the addition of an oxidizing agent to the growth medium. Examples of oxidizing agents include hydrogen peroxide,  $\text{NO}_3^-$ ,  $\text{Fe}^{3+}$ ,  $\text{Ag}^+$ , pyrogallol,  $\text{Cu}^+$  and others. Most preferably, the oxidizing agent will have a redox potential greater than 0.10 and be used at a concentration that avoids toxicity to the growing microorganisms.

The dye-forming substrate is preferably a substituted indigo derivative of the formula:

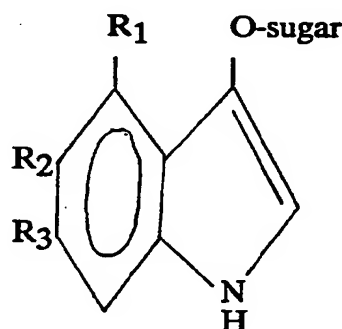

wherein each R substituent is a halo group or hydrogen. The sugar group is any monosaccharide, derivative thereof, or an aldehydo acid. Preferably, the sugar group is a galactoside or a glucoside and the aldehydo acid is glucuronic acid. Most preferably, the sugar or aldehydo acid is attached to oxygen by a beta linkage. Examples of the substituted indigo dyes include 3-indolyl-, 5-bromo-4-chloro-3-indolyl-, 5-bromo-3-indolyl-, 5-bromo-6-chloro-3-indolyl-, and 5-iodo-3-indolyl- derivatives in beta-D linkages with galactose and glucuronic acid. The substrate is colorless, however when hydrolyzed by the appropriate enzyme, the indoxyl or substituted indoxyl portion of the substrate oxidizes to form an intensely colored, insoluble indigo precipitate or indigo derivative as shown, for example, in the Figure. The dimerizing reaction requires an oxidizing environment. Further, non-indigogenic beta-sugar substrates which form colored insoluble products after

enzymatic cleavage could be used for the inventive method and the inventive enzyme indicator device.

In another embodiment, a culture of a particular microorganism or a group of microorganisms can grow in a growth medium or have yielded a presumptive positive result for a particular microorganism or a group of microorganisms. The inventive enzyme indicator device for assaying for a particular enzyme associated with a particular microorganism or group of microorganisms can be added to the growth medium after growth of the sample in the growth medium. For example, a dye-forming substrate can be mixed with a polymer system and then layered on top of a solid support. Alternatively, the dye-forming substrate can be adhered onto the surface of a solid support member. It is preferred that the solid support be a floatable plastic means to remain on the surface of a liquid growth medium. Alternatively, a dye-forming substrate can be mixed in with a polymer material, such as latex, to form a matrix containing the dye-forming substrate throughout. The latex matrix with dye-forming substrate can be dropped into a liquid growth medium or suspended on the surface of the liquid growth medium.

If the dye-forming substrate is not in contact with air, then a developer solution should be added to the growth medium concurrently with the dye-forming substrate or shortly thereafter. The developer solution comprises an oxidizing agent, such as hydrogen peroxide,  $\text{Ag}^+$ ,  $\text{Hg}^{++}$ ,  $\text{I}_3^-$ ,  $\text{IO}_4^-$ , persulfate,  $\text{Pd}^{++}$ , nitrobenzenes, azobenzenes, quinones, 2-chloro-2-nitro-propane,  $\text{K}_3\text{Fe}(\text{CN})_6$ , indolones, and p-hydroxymercaribenzoate. The developer solution preferably further comprises a precipitation enhancer and a cell permeabilizer. The cell permeabilizer facilitates entrance of an indigogenic substrate into the cell to result in a faster turnover of the substrate, and a faster exiting of the indoxyl or indigo from the bacterial cell. Examples of cell permeabilizers include, toluene, anionic and nonionic detergents, and cetyltrimethylammonium bromide. The precipitation enhancer helps to prevent the precipitate from migrating from the solid support and/or the polymer matrix. An example of a precipitation enhancer is an alkaline detergent mixture, such as sodium dodecylsulfate at a 0.1% to 1.0% concentration and a pH of 8.0 to 10.0. The precipitation enhancer is preferably used in a growth medium which shows turbidity or potential quenching of the color reaction.

The enzyme indicator device may be made by dissolving a dye-forming substrate, such as X-Gal, in a solvent such as dimethylformamide (DMF), methyl cellosolve (2-methoxyethanol), or DMSO and mixing with a liquid polymer

system, such as latex. Other adhesive-type polymers that can be dried into a solid or semisolid state include, for example, ethylene vinyl acetate (EVA), polyvinyl alcohol (PVA), and polyvinyl pyrrolidone (PVP). EVA is a polymer soluble in methylene chloride. Additionally, polymer systems that can be liquefied in

5 solvents, such as methylene chloride are useful to liquefy the adhesive-type polymer for mixing with the dye-forming substrate. The mixture of dye-forming substrate and polymer is dried to form a solid or pliable solid. Preferably, the liquid dye-forming substrate and polymer, in a more viscous state, is layered onto the surface of a solid support. Most preferably, a solid support is a porous plastic material that can float on the surface of a liquid.

Alternatively, the dye-forming substrate in solvent can be added directly to the surface of a solid support member and allowed to dry (solvent evaporated). Additionally, the dye-forming substrate can be precipitated out of the solvent, such as DMF, by adding a reagent, such as methylene chloride, in which the dye-forming substrate is insoluble. The precipitate is then layered directly onto a solid support member.

The concentration of dye-forming substrate is sufficient to produce a visible response given the surface area and volume of the reaction surface. For example, the concentration X-Gal on a circular solid support member of polyethylene with a 70 micron pore size of dimensions of one-half inch diameter and one-sixteenth thickness ranges from about 50  $\mu$ g to about 1 mg.

The solid support may float on the surface of a liquid, may reside below the surface of the liquid, or be mechanically placed near the surface of a liquid growth medium by physically attaching it to the culture container or culture container closure. Placing the solid support at or near the liquid-air interface offers close proximity to oxygen in the air to facilitate the oxidation of indoxyl to indigo, thereby improving the speed and sensitivity of the assay. Further, if the solid support and ultimate color localization is located at the top of the culture medium, the reaction is not obscured by the culture broth and it becomes easier to read. If the solid support resides below the surface of the liquid, it is preferred that a developer solution be added to facilitate oxidation of the dye-forming substrate.

When testing for coliform bacteria, common components in the medium, such as in LST broth, can obscure the color reaction from a soluble substrate such as ONPG. Thus, ONPG can only be tested in clear media so that weak reactions can be read. Here, the incorporation of the dye-forming substrate onto a solid support of literally any size can conserve reagent use independent of



sample size. The same solid support and surface area of dye-forming substrate associated with a polymer system can be used for any volume of culture, because the substrate reaction with the bacterial enzyme takes place within a defined surface area rather than throughout the entire volume of the culture. The blue color which is generated is not obscured by ordinary bacteriological media components and therefore an assay with the dye-forming substrate in a polymer system can be performed in a complex medium such as LST broth. The insoluble blue precipitate that is formed from hydrolysis of the dye-forming substrate in the presence of an oxidizing environment remains at the surface of the solid support, thus concentrating the reaction product into a comparatively small, more readily-observed space. Localization of the precipitate makes the test easier to read and increases its sensitivity. The localization of the enzyme assay to a solid support also results in an economy of reagent use.

In another aspect, the inventive method and enzyme indicator device can be used in a quantitative assay to estimate the number of coliform bacteria in a food or water sample using the most probable number (MPN) method. An MPN coliform assay is performed by incorporating a dye-forming substrate in a polymer system, such as latex into a standard coliform medium, such as LST broth. Preparation of the sample is done by standard techniques. After overnight incubation, the tubes are read for color production in the polymer system and the positive tubes are used to calculate a confirmed MPN per gram or MPN per milliliter value, according to standard MPN techniques. See Compendium of Methods for the Microbiological Examination of Foods (2nd Ed., American Public Health Association, 1984), Official Methods of Analysis (14th Ed., Association of Official Analytical Chemists, 1984), and Standard Methods for the Examination of Water and Wastewater (17th Ed., American Public Health Association, 1981).

The inventive method and enzyme indicator device are used, for example, for qualitative testing for coliform bacteria. An indigogenic galactoside substrate is incorporated into the medium used for the PA test for coliforms in drinking water. The water sample is added to a medium containing the dye-forming substrate in a polymer system, preferably layered onto a floatable plastic solid support. Precipitated indigo, on the solid support member, within the culture vessel after a suitable incubation period signifies the presence of coliforms in the sample. This procedure is also used for qualitative testing of environmental samples for the presence of coliforms. For example, a surface or drain could be assayed for the effectiveness of sanitation by swabbing with a cotton-tipped

applicator and then placing the applicator in a coliform medium containing the indigogenic galactoside dye-forming substrate incorporated into a polymer system, preferably on a solid support member. A precipitated indigo color within the culture vessel, or more specifically on the surface of the solid support after a suitable incubation period, such as 1 to 24 hours at a permissive temperature for bacterial growth would signify the presence of coliforms in the sample.

In addition to the use of a dye-forming substrate as described above, the present invention provides for the simultaneous confirmatory determination of both E. coli and coliform bacteria without the need for additional culture transfers, reagent preparation or diagnostic manipulating. More specifically, a disc as described above containing dye-forming and fluorogenic-generating substrates is added to a container of microbiological culture media to be tested for the presence of E. coli and/or the coliform group of bacteria. The container with the medium and reaction disc is incubated at a permissive temperature for growth of the microorganisms of interest. Following an incubation step, the container is observed for the presence of a distinguishing color on the reaction disc or in the surrounding medium. The container is then exposed to long wave ultraviolet light (366 nm) and observed for fluorescence in the surrounding medium. The observance of color on the disc indicates a confirmed reaction for the coliform group of organisms. The presence of fluorescence indicates a confirmed reaction for the presence of the specific microorganism, E. coli.

The following examples are offered by way of illustration and not by way of limitation.

#### EXAMPLE 1

This example illustrates a method for detecting coliform bacteria in drinking water by the PA (presence-absence) technique, as described in Standard Methods for the Examination of Water and Wastewater, (17th Ed., American Public Health Association, 1981). Briefly, 100 ml of water are aseptically added to 50 ml of triple strength Clark's medium in a 250 ml bottle. Clark's medium had been modified to exclude a pH indicator. After inoculation, a plastic support containing an indigogenic galactoside substrate in a latex polymer system was added to the culture vessel so that it floated at the liquid surface. The indigogenic galactoside substrate in a latex polymer system on the plastic support was made as described herein. The mixture was incubated at 35°C and the plastic support was inspected for indigo formation at 24 hours. The indigo-positive reaction signified the confirmed presence of coliforms in the water sample.

### EXAMPLE 2

This example illustrates a method for detecting coliform bacteria and E. coli in an environmental sample using the inventive method. The environmental sample was taken in an effort to assess the effectiveness of sanitation or as part of a routine environmental screening program. A cotton-tipped applicator was wetted with sterile diluent and used to swab a surface. The swab was placed in a bacteriological lauryl sulfate tryptose broth growth medium. A buoyant circular plastic support layered with an indigogenic galactoside substrate in a polymer system on one-half of the circle and an indigogenic glucuronide substrate on the other half of the circle was added after inoculation. The mixture was incubated for 35°C for 24 hours. An indigo-positive reaction on one half of the plastic support signified the presence of coliform bacteria in the environmental sample. An indigo-positive reaction on both halves of the plastic support signified the presence of E. coli in the sample.

### EXAMPLE 3

This example illustrates a method for estimating the level of coliforms and E. coli in a food or water sample by the Most Probable Number (MPN) analysis using the inventive method. The MPN technique is described in Compendium of Methods for the Microbiological Examination of Foods (2nd Ed., American Public Health Association, 1984), Official Methods of Analysis (14th Ed., Association of Official Analytical Chemists, 1984), and Standard Methods for the Examination of Water and Wastewater (17th Ed., American Public Health Association, 1981).

Water samples (10 ml) were aseptically added to 10 ml of double-strength lauryl sulfate tryptose (LST) broth. A buoyant circular plastic solid support containing an indigogenic galactoside substrate in a polymer system on one-half of the circular plastic and an indigogenic glucuronide substrate in a polymer system on the other half was added to each test tube after inoculation. The mixture was incubated at 35° C for 24 hours. An indigo-positive reaction on one half of the plastic support signifies the confirmed presence of coliform bacteria in the water sample. An indigo-positive reaction on both halves of the plastic support confirms the presence of E. coli in the sample. The number of indigo-positive cultures was used to estimate the number of coliform bacteria and E. coli present within the original water sample through the use of a standard statistical probability table (MPN Table).

The food samples were diluted in serial 10-fold fashion prior to the MPN assay. One part (by mass) of solid food sample was added to 9 parts (by volume) of sterile diluent Butterfield's buffered phosphate buffer. The mixture was blended at high speed for two minutes. Two additional 10-fold dilutions were made from the 1:10 diluted food sample. One ml of each 10-fold dilution was added to each of three 10 ml test tubes of LST broth. A buoyant circular plastic support containing an indigogenic galactoside substrate in a polymer system on one half and an indigogenic glucuronide substrate in a polymer system on the other half was added to each test tube after inoculation. The mixture was incubated at 35°C for 24 hours. An indigo-positive reaction on one half of the plastic support confirmed the presence of coliform bacteria in the sample dilution. An indigo-positive reaction on both halves of the plastic support confirmed the presence of E. coli in the sample dilution. The number of indigo-positive cultures was used to estimate the number of coliform bacteria and E. coli present within the original food sample through the use of standard statistical probability tables (MPN Table).

#### EXAMPLE 4

This example illustrates a method for the confirmation of presumptively positive coliform cultures for total coliforms and for E. coli using an indigogenic substrate mixture. Cultures which were presumptively gas-positive for coliforms were confirmed for coliform and/or E. coli presence according to the following procedure. A buoyant circular solid support containing an indigogenic galactoside substrate on one half and an indigogenic glucuronide substrate on the other half was added to the culture. An oxidizing agent (hydrogen peroxide at 10% concentration) was added to the culture. The presumptive portion of the coliform MPN analysis was performed in LST broth on a food or water sample using the procedure described in Example 3. Cultures which were gas-positive after 24 to 48 hours of incubation were confirmed with the buoyant, circular solid support described herein along with the addition of a developer comprising the oxidizing agent. The culture was incubated at 35°C for approximately one hour and the plastic support was read for indigo formation. Cultures which were indigo-positive on one half of the circular solid support were considered to be confirmed for the presence of coliform bacteria. Cultures which were indigo-positive on both halves of the circular solid support were confirmed for the presence of E. coli. Negative cultures were reincubated for an additional hour

and reread. Cultures which were still negative for indigo at two hours remained unconfirmed for the presence of coliforms or E. coli.

#### EXAMPLE 5

5           This example illustrates a method for detecting coliform bacteria in drinking water by the PA (presence-absence) technique, as described in Standard Methods for the Examination of Water and Wastewater, (17th Ed., American Public Health Association, 1981). Briefly, 100 ml of water are aseptically added to 50 ml of triple strength Clark's medium in a 250 ml bottle. Clark's medium had  
10       been modified to exclude a pH indicator. After inoculation, a plastic support containing an indigogenic galactoside substrate and a fluorogenic-generating substrate in a latex polymer system was added to the culture vessel so that it floated at the liquid surface. The indigogenic galactoside substrate in a latex polymer system on the plastic support was made as described herein. The mixture  
15       was incubated at 35°C and the plastic support was inspected for indigo formation at 24 hours. The indigo-positive reaction signified the confirmed presence of coliforms in the water sample.


#### EXAMPLE 6

20           This example illustrates a method for detecting coliform bacteria and E. coli in an environmental sample using the inventive method. The environmental sample was taken in an effort to assess the effectiveness of sanitation or as part of a routine environmental screening program. A cotton-tipped applicator was wetted with sterile diluent and used to swab a  
25       surface. The swab was placed in a bacteriological lauryl sulfate tryptose broth growth medium. A buoyant circular plastic support layered with an indigogenic galactoside substrate in a polymer system and a fluorogenic glucuronide substrate was added after inoculation. The mixture was incubated for 35°C for 24 hours. An indigo-positive reaction on the plastic support signified the presence of  
30       coliform bacteria in the environmental sample. A fluorescent-positive reaction in the container when received under long wave ultraviolet light (366 nm) of the plastic support signified the presence of E. coli in the sample.

          Although the foregoing invention has been described, in part, by way of illustration and example for the purposes of clarity and understanding, it  
35       will be apparent that certain changes or modifications may be practiced without deviating from the spirit and scope of the invention.

Claims**We claim:**

1. A method for assaying for an enzyme associated with a particular microorganism or a group of microorganisms comprising:
  - adding a sample to a liquid growth medium;
  - adding a dye-forming substrate on a defined reaction surface wherein the dye-forming substrate forms a precipitate when cleaved by the enzyme;
  - incubating the sample, growth medium and dye-forming substrate at a permissive temperature for the growth of the particular microorganism or group of microorganisms; and
  - detecting the presence of the particular microorganism or group of microorganisms by a colored precipitate on the defined reaction surface.
2. The method of claim 1 wherein the dye-forming substrate forms an insoluble precipitate upon enzymatic cleavage.
3. The method of claim 2 wherein the dye-forming substrate is a substituted indigo derivative of the formula:

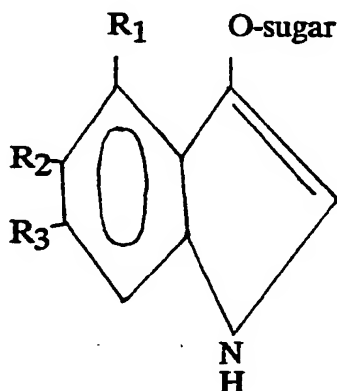

wherein each R substituent is a halo group or hydrogen.

4. The method of claim 3 wherein the halo group is selected from the group consisting of chloro, bromo, and iodo.

5. The method of claim 3 wherein the sugar is a monosaccharide or derivative thereof.

6. The method of claim 5 wherein the sugar is selected from the group consisting of glucuronide, glucoside, and galactoside.

7. The method of claim 1 wherein the dye-forming substrate further comprises a polymer system, wherein the dye-forming substrate is incorporated in the polymer system to form the reaction surface.

8. The method of claim 7 further comprising a solid support member, wherein the dye-forming substrate in the polymer system is layered on a surface of the solid support member, and wherein the surface of the solid support member becomes the reaction surface.

9. The method of claim 8 wherein the solid support is a plastic member that floats on the surface of a liquid and allows the dye-forming substrate in a polymer system to be in contact with the air above the liquid.

10. The method of claim 8 wherein the solid support is a side of a culture vessel that is in contact with both the liquid and the air above the liquid.

11. The method of claim 8 wherein the solid support is physically suspended to reside at the liquid-air interface.

12. The method of claim 7 wherein the polymer system is selected from the group consisting of latex, ethylene vinyl acetate, polyvinyl alcohol, polyvinyl pyrrolidone and combinations thereof.

13. The method of claim 12 wherein the polymer system is a latex.

14. The method of claim 1 wherein the particular microorganism is E. coli and the group of microorganisms is the coliform group of bacteria.

15. The method of claim 1 wherein the sample is a water sample, a drinking water sample, waste water sample, a food sample, or an environmental sample.

16. The method of claim 1 further comprising adding a developer solution during the incubating step, wherein the developer solution comprises an oxidizing agent.

17. The method of claim 16 wherein the oxidizing agent is selected from the group consisting of hydrogen peroxide,  $\text{Ag}^+$ ,  $\text{Hg}^{++}$ ,  $\text{I}_3$ ,  $\text{IO}_4$ ,  $\text{Pd}^{++}$ , persulfate, p-hydroxymercuribenzoate, and combinations thereof.

18. The method of claim 16 further comprising adding a precipitation enhancer and a cell permeabilizer to the developer solution.

19. The method of claim 18 wherein the cell permeabilizer is selected from the group consisting of toluene, detergent, cetyltrimethyl ammonium bromide and combinations thereof.

20. The method of claim 18 wherein the precipitation enhancer is selected from the group consisting of alkaline detergents.

21. A method for assaying for an enzyme associated with a particular microorganism or group of microorganisms comprising:

adding a sample to a liquid growth medium;

incubating the sample and growth medium at a permissive temperature to allow growth of the particular microorganism or group of microorganisms until there is evidence of microorganism growth;

adding an enzyme indicator device, wherein the enzyme indicator device comprises a dye-forming substrate on a defined reaction surface, and wherein the dye-forming substrate is cleaved by the enzyme to form a colored precipitate;


incubating the sample growth medium and enzyme indicator device at a permissive temperature from about 1 hour to about 48 hours; and

visually detecting the presence of the particular microorganism or group of microorganisms by a color change on the enzyme indicator device.

22. The method of claim 21 wherein the dye-forming substrate is a substituted indigo derivative of the formula:



-23-

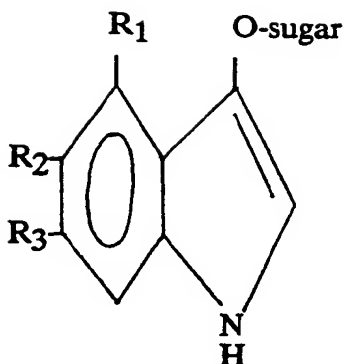

wherein each R substituent is a halo group or hydrogen.

23. The method of claim 22 wherein the halo group is selected from the group consisting of chloro, bromo, and iodo.

24. The method of claim 21 wherein the enzyme indicator device further comprises a polymer system, wherein the dye-forming substrate is incorporated in the polymer system to form the reaction surface.

25. The method of claim 24 wherein the polymer system is selected from the group consisting of latex, ethylene vinyl acetate, polyvinyl alcohol, polyvinyl pyrrolidone and combinations thereof.

26. The method of claim 24 wherein the polymer system is a latex.

27. The method of claim 24 wherein the enzyme indicator device is formed by dissolving the dye-forming substrate in dimethylformamide, mixing the dye-forming substrate in liquid polymer system to form a mixture, and drying the layered solid support.

28. The method of claim 24 wherein the enzyme indicator device further comprises a solid support, wherein the dye-forming substrate incorporated in the polymer system is layered on a surface of the solid support to form the reaction surface.

29. The method of claim 28 wherein the solid support is clear, turbid, or white colored.

30. The method of claim 21 wherein the defined reaction surface comprises a solid support member.

31. The method of claim 21 wherein the particular microorganism is E. coli and the group of microorganisms is the coliform group of bacteria.

32. The method of claim 21 wherein the sample is a water sample, a waste water sample, a drinking water sample, a food sample, or an environmental sample.

33. The method of claim 21 wherein the enzyme indicator device is added to the top surface of the growth medium.

34. The method of claim 21 further comprising adding a developer solution together with the enzyme indicator device, wherein the developer solution comprises an oxidizing agent.

35. The method of claim 34 wherein the oxidizing agent is selected from the group consisting of hydrogen peroxide,  $\text{Ag}^+$ ,  $\text{Hg}^{++}$ ,  $\text{I}_3$ ,  $\text{IO}_4$ ,  $\text{Pd}^{++}$ , persulfate, p-hydroxymercuribenzoate, and combinations thereof.

36. The method of claim 34 further comprising adding a precipitation enhancer and a cell permeabilizer to the developer solution.

37. The method of claim 36 wherein the cell permeabilizer is selected from the group consisting of toluene, detergent, cetyltrimethyl ammonium bromide and combinations thereof.

38. The method of claim 36 wherein the precipitation enhancer is selected from the group consisting of alkaline detergents.

39. An enzyme indicator device for assaying for a particular microorganism or group of microorganisms associated with an enzyme, comprising a dye-forming substrate on a defined reaction surface, wherein the dye-forming substrate is cleaved by the enzyme to form a colored precipitate on the reaction surface.


40. The enzyme indicator device of claim 39 further comprising a solid support member, wherein the solid support member defines the reaction surface.

41. The enzyme indicator device of claim 40 wherein the solid support floats on the surface of a liquid.

42. The enzyme indicator device of claim 40 wherein the solid support is suspended at the liquid-air interface.

43. The enzyme indicator device of claim 40 wherein the reaction surface is incorporated on a culture vessel.

44. The enzyme indicator device of claim 39 wherein the dye-forming substrate is a substituted indigo derivative of the formula:

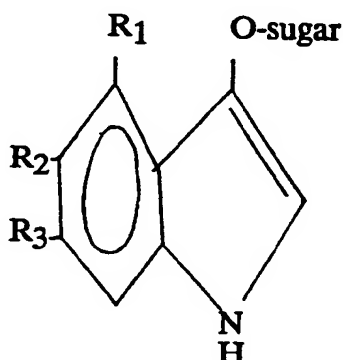

wherein each R substituent is a halo group or hydrogen.

45. The enzyme indicator device of claim 44 wherein the halo group is selected from the group consisting of chloro, bromo, and iodo.

46. The enzyme indicator device of claim 39 further comprising a polymer system, wherein the dye-forming substrate is incorporated in the polymer system to form the reaction surface.


47. The enzyme indicator device of claim 46 wherein the polymer system is selected from the group consisting of latex, ethylene vinyl acetate, polyvinyl alcohol, polyvinyl pyrrolidone and combinations thereof.

48. The enzyme indicator device of claim 47 wherein the polymer system is a latex.

49. The enzyme indicator device of claim 39 wherein the particular microorganism is E. coli and the group of microorganisms is the coliform group of bacteria.

50. An enzyme indicator device for assaying for a particular microorganism or group of microorganisms comprising a dye-forming substrate dried on the surface of a solid support member.

51. The enzyme indicator device of claim 50 wherein the dye-forming substrate is a substituted indigo derivative of the formula:

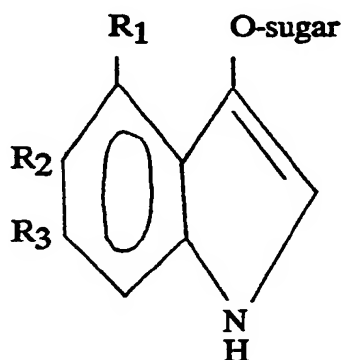

wherein each R substituent is a halo group or hydrogen.

52. The enzyme indicator device of claim 51 wherein the halo group is selected from the group consisting of chloro, bromo, and iodo.

53. The enzyme indicator device of claim 50 wherein the particular microorganism is E. coli and the group of microorganisms is the coliform group of bacteria.

54. A method for assaying for an enzyme associated with a particular microorganism and a group of microorganisms comprising:  
adding a sample to a liquid growth medium;


adding a dye-forming and/or fluorogenic-generating substrate on a defined reaction surface wherein the substrate forms a precipitate or a fluorescent by-product when cleaved by the enzyme;

incubating the sample, growth medium and substrate at a permissive temperature for the growth of the particular microorganism or group of microorganisms; and

detecting the presence of the particular microorganism and group of microorganisms by a colored precipitate and fluorescence on the defined reaction surface or in the surrounding medium.

55. The method of claim 54 wherein the fluorogenic-generating substrate is 4-methylumbelliferyl- $\beta$ -D-glucuronide.

1/1

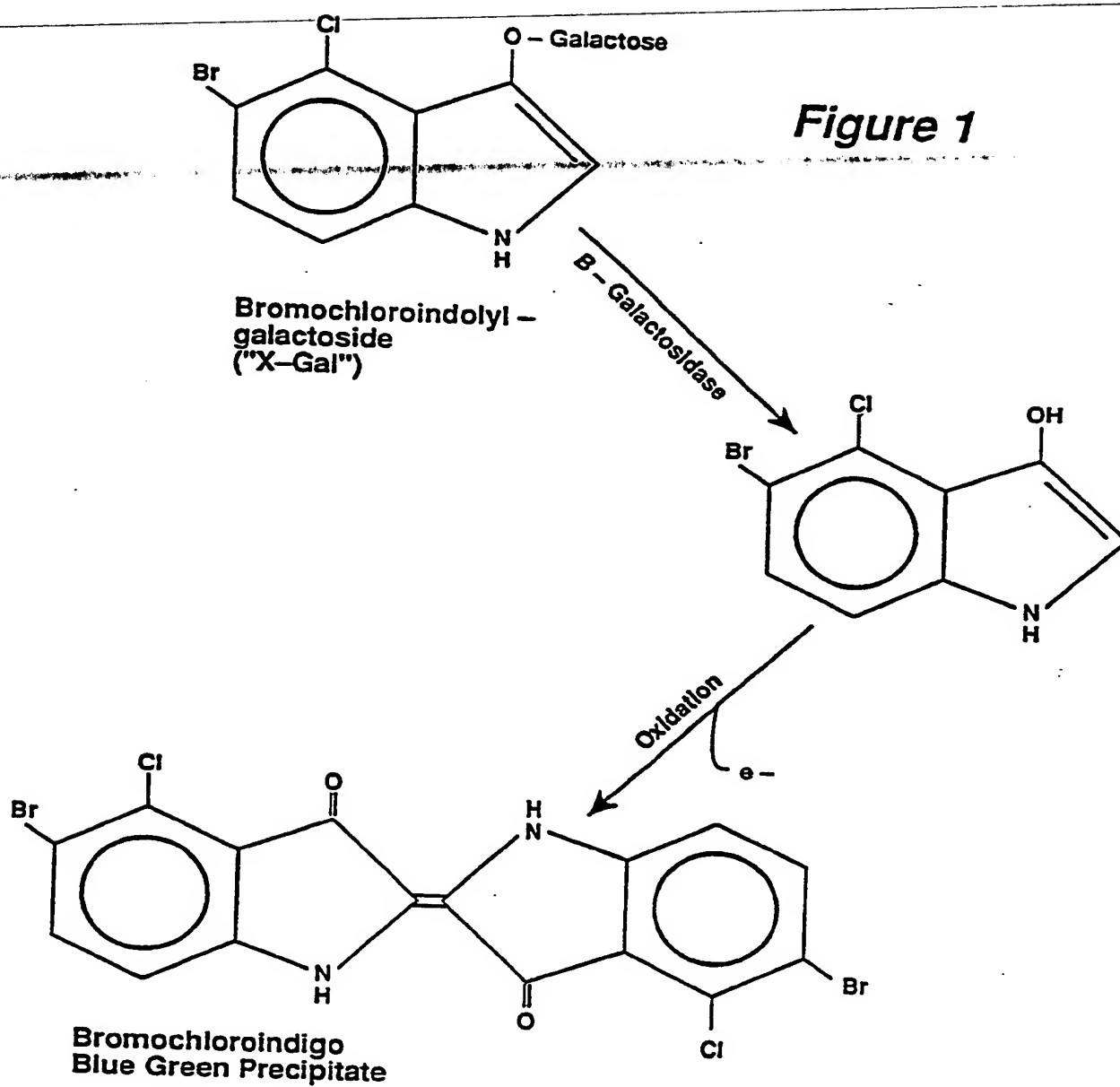

# INTERNATIONAL SEARCH REPORT

International Application No PCT/US 90/02346

## I. CLASSIFICATION OF SUBJECT MATTER (if several classification symbols apply, indicate all) \*

According to International Patent Classification (IPC) or to both National Classification and IPC

IPC<sup>5</sup>: C 12 Q 1/04, C 12 Q 1/34

## II. FIELDS SEARCHED

Minimum Documentation Searched ?

Classification System

Classification Symbols

IPC<sup>5</sup>

C 12 Q

Documentation Searched other than Minimum Documentation  
to the Extent that such Documents are Included in the Fields Searched \*

## III. DOCUMENTS CONSIDERED TO BE RELEVANT \*

Category *	Citation of Document, <sup>11</sup> with indication, where appropriate, of the relevant passages <sup>12</sup>	Relevant to Claim No. <sup>13</sup>
X	EP, A, 0122028 (E-Y LABORATORIES, INC.) 17 October 1984 see the whole document  --	1,3-8,14, 21-24, 28-33,39, 40,44-46, 49-55
Y	EP, A, 0282733 (HER MAJESTY THE QUEEN) 21 September 1988 see the whole document  --	1-55
Y	EP, A, 0128521 (EASTMAN KODAK COMPANY) 19 December 1984 see abstract; pages 1-3; pages 14-15  --	1-55
A	EP, A, 0238243 (COGENT LIMITED) 23 September 1987  -----	

\* Special categories of cited documents: <sup>10</sup>

"A" document defining the general state of the art which is not considered to be of particular relevance

"E" earlier document but published on or after the international filing date

"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)

"O" document referring to an oral disclosure, use, exhibition or other means

"P" document published prior to the international filing date but later than the priority date claimed

"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

"X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step

"Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.

"A" document member of the same patent family

## IV. CERTIFICATION

Date of the Actual Completion of the International Search

2nd August 1990

Date of Mailing of this International Search Report

12. 09. 90

International Searching Authority

EUROPEAN PATENT OFFICE

Signature of Authorized Officer

M. SOTELO

**ANNEX TO THE INTERNATIONAL SEARCH REPORT  
ON INTERNATIONAL PATENT APPLICATION NO.**

US 9002346

SA 37033

This annex lists the patent family members relating to the patent documents cited in the above-mentioned international search report. The members are as contained in the European Patent Office EDP file on 05/09/90. The European Patent Office is in no way liable for these particulars which are merely given for the purpose of information.

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
EP-A- 0122028	17-10-84	JP-A- 59220199	11-12-84
EP-A- 0282733	21-09-88	JP-A- 64002596	06-01-89
		US-A- 4923804	08-05-90
EP-A- 0128521	19-12-84	US-A- 4556636	03-12-85
		CA-A- 1211349	16-09-86
		JP-A- 60006199	12-01-85
EP-A- 0238243	23-09-87	WO-A- 8705632	24-09-87

EPO FORM P0479

For more details about this annex : see Official Journal of the European Patent Office, No. 12/82